

Genuair


TECHNIQUE

1. Remove the protective cap by **lightly squeezing the arrows** marked on each side and pulling outwards.
2. Hold the inhaler **horizontally** with the mouthpiece towards you and the colored button facing **straight up**.
3. **Press** the colored button all the way down and **release** it. The **green** control window confirms that your medicine is ready for inhalation.
4. Breathe out away from the inhaler.
5. Place your lips tightly around the mouthpiece and inhale **strongly** and **deeply** through your mouth. The inhaler signals that you inhaled correctly by emitting a **“click”** sound. Keep breathing in even after you have heard the **“click”** to be sure you get the full dose.
6. Remove the inhaler from your mouth and hold your breath for 5-10 seconds or as long as is comfortable, then breathe out slowly.
7. Make sure the control window has turned to **red**. This confirms that you have inhaled your full dose correctly. If the colored control window is still green, repeat inhaling strongly and deeply through the mouthpiece. If the window does not change to red, you may have forgotten to release the green button before inhaling, or may have not inhaled correctly. If that happens, try again.
8. Replace the protective cap by pressing it back onto the mouthpiece.

www.livingwellwithcopd.com

Genuair

TECHNIQUE -continued

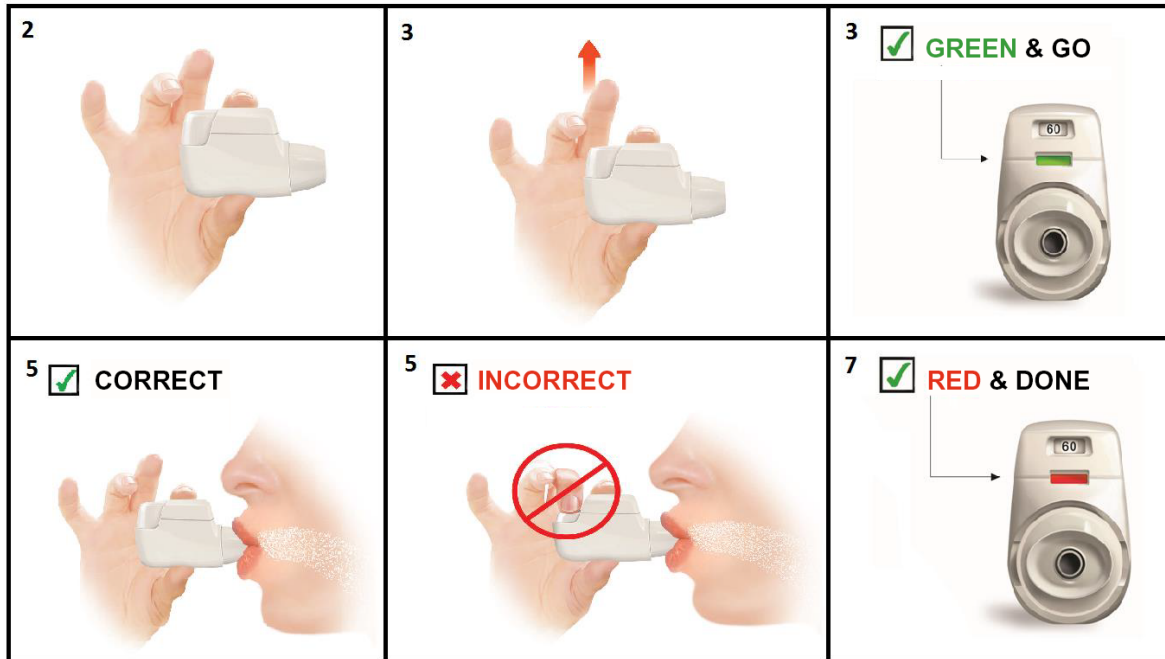

Maintenance and Cleaning

- When a red striped band appears in the dose indicator this means you are nearing your last dose and you should obtain a new inhaler.
- When 0 (zero) appears in the middle of the dose indicator, you should continue using any doses remaining in the inhaler.
- When the last dose has been prepared for inhalation, the green button will not return to its full upper position, but will be locked in a middle position. Even though the green button is locked, your last dose may still be inhaled. After that, the inhaler cannot be used again and you should start using a new inhaler.
- You do not need to clean your inhaler. Never use water to clean the inhaler, as this may damage your medicine.

www.livingwellwithcopd.com